



Finite element analysis regarding the propagation of transient vibrations from handheld equipment

The sample measured is a replica of the skin on a fingertip. In this case, the difficulty of holding one's finger absolutely still during the measurement necessitates using the indirect method of measuring on replicas. The objective was to obtain information for simulation of microgeometry of skin to be able to create a simplified virtual model. The model was going to be used in finite element analysis (FEA) regarding the propagation of transient vibrations from handheld equipment, such as impact wrenches and chisel machines, through the skin to the nerves where damage occurs. These vibrations can cause a medical condition called "white fingers" which is a common work-related health problem.


Before this government-sponsored project, the microgeometry of skin had not yet been characterized for this specific purpose. This project's main objective was getting detailed information regarding skin's microgeometry over a large enough area to be able to draw general conclusions for use in further investigations.




[Swerea IVF](#) is a research institute that offers advanced research and consulting services to put new technology and new methods into practice with a focus on product, process and production development. Swerea IVF also offers in-depth expertise in relation to materials' properties and applications for ceramic, polymer and textile materials. Swerea IVF is part of the Swerea group, a Swedish research group for industrial renewal and sustainable development. They produce, develop, and disseminate research results in the areas of materials development, production development and product development. Project current being undertaken by Johan Berglund and Hans Lindell.

swerea
swedish research


■ Measurements

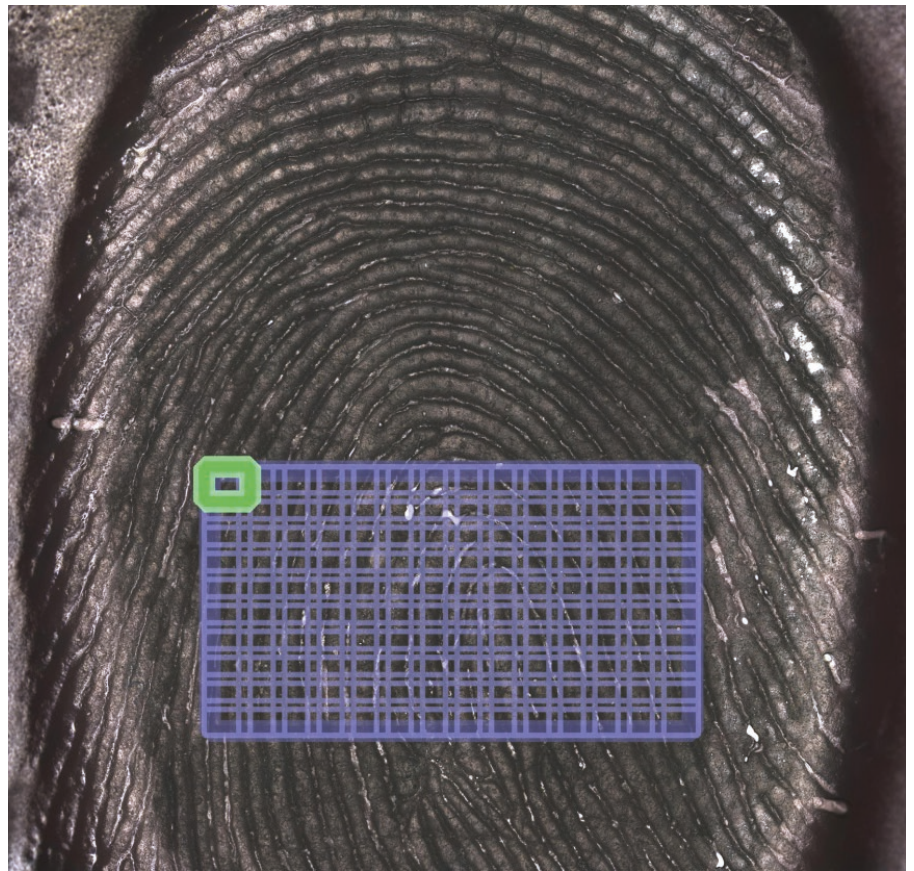

Overview of skin with the measurement area depicted inside the rectangular blue grid.

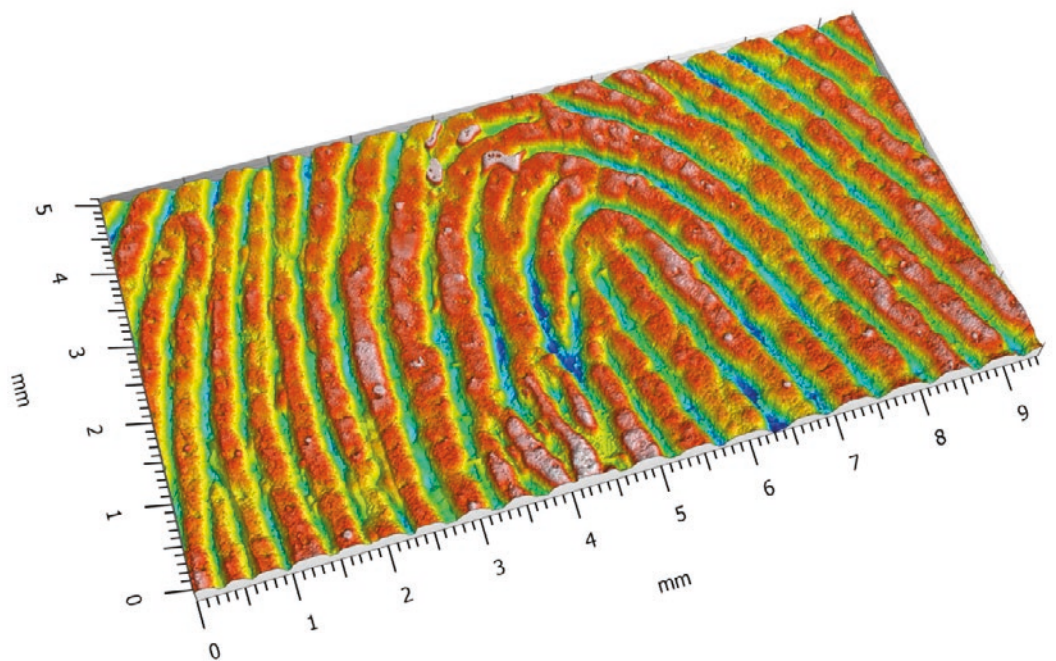

5x9 mm Extended measurement.

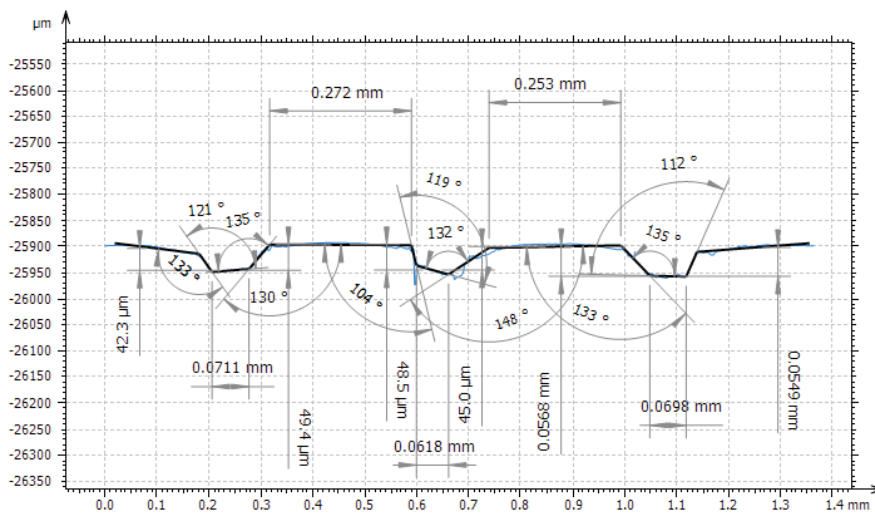

Example of an evaluation of the sample's geometry using SensoMAP advanced contour module.

Conclusions

Using the confocal technique, the 20X brightfield objective lens, and extended measurements, we were able to measure the skin geometry with a high enough resolution and sufficient reliability to draw general conclusions regarding the size, height, and distribution of the ridges and lines of the skin. With this data, we were able to make a simplified model of the skin for successful use in FEA.

For more information, check the [Attenuation of hand-held machine vibrations. Application of non-linear tuned vibration absorber.](#)



SENSOFAR is a leading-edge technology company that has the highest quality standards within the field of surface metrology

Sensofar Metrology provides high-accuracy optical profilers based on confocal, interferometry and focus variation techniques, from standard setups for R&D and quality inspection laboratories to complete non-contact metrology solutions for in-line production processes. The Sensofar Group has its headquarters in Barcelona, also known as a technology and innovation hub in Europe. The Group is represented in over 30 countries through a global network of partners and has its own offices in Asia, Germany and the United States.

HEADQUARTERS

SENSOFAR METROLOGY | BARCELONA (Spain) | T. +34 93 700 14 92 | info@sensofar.com

SALES OFFICES

SENSOFAR ASIA | SHANGHAI (China) | T. +86 021 51602735 | info.asia@sensofar.com

SENSOFAR GERMANY | MUNICH (Germany) | T. +49 151 14304168 | info.germany@sensofar.com

SENSOFAR USA | NEWINGTON (USA) | T. +1 617 678 4185 | info.usa@sensofar.com

sensofar.com